



**Quality You Can Feel.
Safety You Can Trust**

www.larsmedicare.com 

We are a world-class manufacturer of comprehensive range of premium medical disposable devices, including IV Cannula, Safety IV Cannula, Three Way Stop Cock With/Without Extension Tube, IV Flow Regulator, Infusion Set, Blood Transfusion Set, Measured Volume Set, Syringes, Latex Gloves, Respiratory and much more.

With state-of-the-art manufacturing in Haryana & Gujarat and head office in New Delhi & Americas Headquarters in Canada, we proudly deliver ISO-certified, CE-marked products to over 100 countries worldwide.

Lars Medicare – Precision, Safety & Care in Every Product.



Our Mission

To manufacture world-class medical disposables ensuring safety, reliability, and accessibility for global healthcare needs.



Our Vision

To be a trusted global leader in medical disposables, empowering healthier lives through innovation and care.



Our Values

Integrity, innovation, compassion, quality, and commitment drive our pursuit of excellence in healthcare solutions.

At Lars Medicare, our journey of over two decades has been guided by a single vision to deliver reliable, innovative, and world-class medical disposables that empower healthcare professionals and improve patient outcomes worldwide.

From humble beginnings to becoming a globally trusted name, our success is rooted in a commitment to quality, safety, and innovation. Every product we design and manufacture reflects our dedication to meeting the highest international standards, ensuring that caregivers can work with confidence and patients can receive the care they deserve.

We believe that healthcare is not just about technology but about trust. This is why we continue to invest in advanced manufacturing, research, and sustainable practices, ensuring that we remain at the forefront of medical innovation while being mindful of our responsibility towards people and the planet.

As we expand our global footprint, our focus remains clear: to build long-term partnerships, provide exceptional value, and contribute meaningfully to the future of healthcare.

On behalf of the entire Lars Medicare family, I thank our partners, clients, and team members for their trust and continued support. Together, we will keep striving to create a healthier and safer tomorrow.

Sanjay Goel
Managing Director

Our Certification




IV Cannula

Provein™

With Wings and Injection Valve

Available Size: 14, 16, 17, 18, 20, 22, 24, 26G

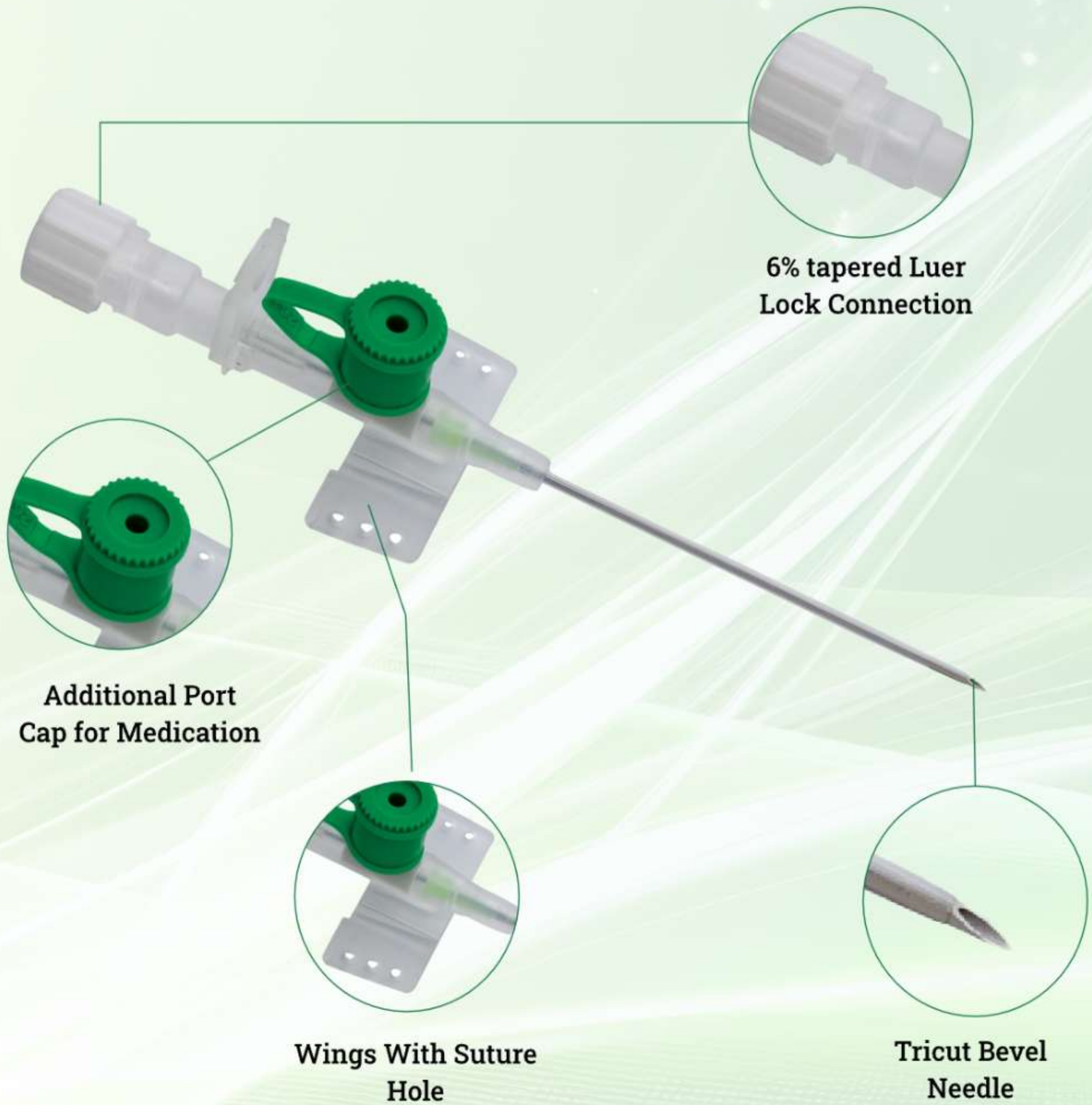

Provein α TM

Without Wings and Injection Valve.
Available Size: 14, 16, 17, 18, 20, 22, 24, 26G



Provein β TM

With Wings and without Injection Valve.
Available Size: 14, 16, 17, 18, 20, 22, 24, 26G

PropenTM

Pen Type IV Cannula without Wings and Injection Valve.

With long and broad needle hub for a perfect grip.

Available Size: 18, 20, 22, 24G



AVAILABLE IN PTFE/PUR/FEP CATHETER

Technical Specification

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
14	Orange	1.6-2.1	45	270
16	Grey	1.3-1.7	45	180
17	White	1.2-1.5	45	125
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	32	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20
26	Violet	0.4-0.6	19	13

Packaging Details

Product	No. of Pieces/ Carton	Carton Dimension (cm)		
		L	W	H
Provein (14G 16G 17G 18G)	1000	70.00	32.00	18.00
Provein (20G 22G 24G 26G)	1000	63.00	32.00	18.00
Provein α	1000	54.00	29.00	18.00
Provein β	1000	70.00	32.00	18.00
Pro-Pen	500	32.00	29.00	10.00

Safety IV Cannula with Plastic Cage

Provein Safe™

- With Wings and Injection Valve.
- Available Size: 14, 16, 17, 18, 20, 22, 24G



6% Tapered Luer Lock Connection



Flash Back Chamber



Flip Port Cap



Plastic Cage Safety Shield

Provein Safe α TM

Without Wings and Injection Valve.
Available Size: 14, 16, 17, 18, 20, 22, 24G



Provein Safe β TM

With wings and without injection valve.
Available Size: 14, 16, 17, 18, 20, 22, 24G



AVAILABLE IN PTFE/PUR CATHETER

Technical Specification

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
14	Orange	1.6-2.1	45	270
16	Grey	1.3-1.7	45	180
17	White	1.2-1.5	45	125
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	32	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20

Packaging Details


Product	No. of Pieces /Carton	Carton Dimension (in cm)		
		L	W	H
Provein Safe	500	45.00	36.00	18.00
Provein safe α	500	45.00	36.00	18.00
Provein Safe β	500	45.00	36.00	18.00

Safety IV Cannula with Metallic Clip



Safe Inject

Wings and Injection Valve.
Available Size: 18, 20, 22, 24G

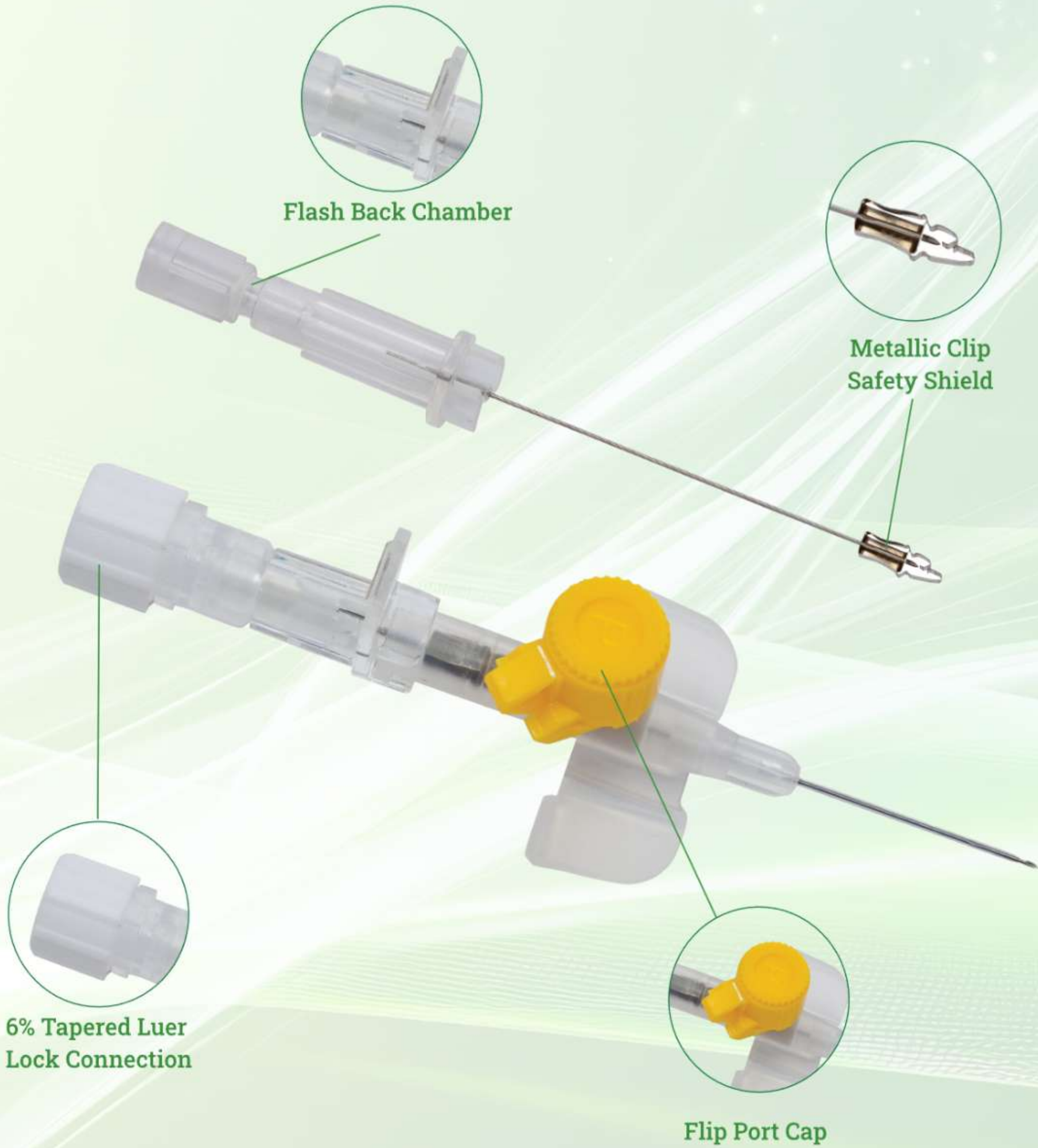

Safe Cath™

Without Wings and Injection Port
Available Size: 18, 20, 22, 24G



Safe Cath Pro™

Without Wings and Injection Port.
With Long Projection Body & Flash Back Chamber.
Available Size: 18, 20, 22, 24G



AVAILABLE IN PTFE/PUR CATHETER

Technical Specification

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water Flow (ml/min)
18	Green	0.9-1.3	45	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Safe Inject	500	45.00	36.00	18.00
Safe Cath	500	45.00	36.00	18.00
Safe Cath Pro	500	45.00	36.00	18.00


Blood Control Technology (BCT)

Blood control technology minimize the risk of blood exposure while doing cannulation as mechanical valve closes when the needle is withdrawn. There is no venous compression required during vein puncture and increases patient and healthcare worker safety.

Provein 2.0™

IV Cannula with Wings and Injection Valve and with **Blood Control Technology**

Available Size: 18, 20, 22, 24G

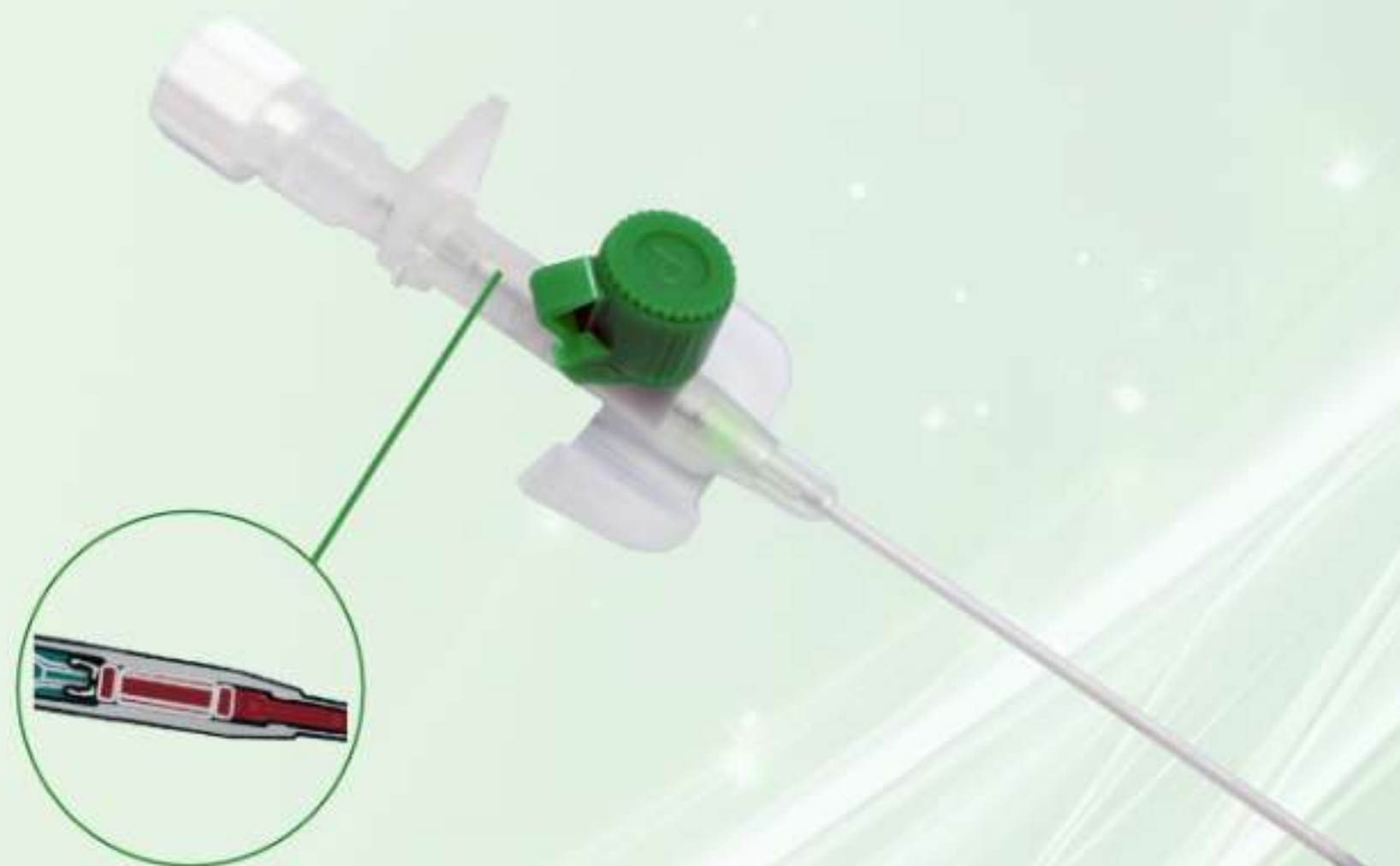

Blood Control Valve

Provein α 3.0™

IV Cannula without Wings and Injection Valve and with **Blood Control Technology**

Available Size: 18, 20, 22, 24G



Blood Control Valve

Technical Specification

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	32	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20

Packaging Details

Product	No. of Pieces/ Carton	Carton Dimension (cm)		
		L	W	H
Provein 2.0	500	45.00	36.00	18.00
Provein α 3.0	500	45.00	36.00	18.00

Provein Safe - BCT

Safety IV Cannula with Wings and Injection Valve and with **Blood Control Technology**

Available Size: 18, 20, 22, 24G



Safe Cath - BCT

Safety IV Cannula without Wings and Injection Valve and with **Blood Control Technology**

Available Size: 18, 20, 22, 24G



AVAILABLE IN PTFE/PUR CATHETER

Technical Specification

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water Flow (ml/min)
18	Green	0.9-1.3	45	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20

Packaging Details


Product	No. of Pieces /Carton	Carton Dimension (in cm)		
		L	W	H
Provein Safe - BCT	500	45.00	36.00	18.00
Safe Cath - BCT	500	45.00	36.00	18.00

Closed Loop IV Cannula with BCT

Gen-X™

Closed Loop IV Cannula with Blood Control Technology

- All-in-one solution for peripheral venous access by combining the catheter, extension set, stabilization device & needle-free connector in single unit.
- Reduces infection risk & related complications.
- Blood Control Technology minimize the risk of blood exposure.
- Facilitate triple infusion (3 ports device).
- Safety mechanism prevents needle stick injury thus provide safety to patient and healthcare worker.
- Polyurethane catheter material ensures better bio-compatibility & longer indwelling period.
- Flashback window for quick visualization of blood.

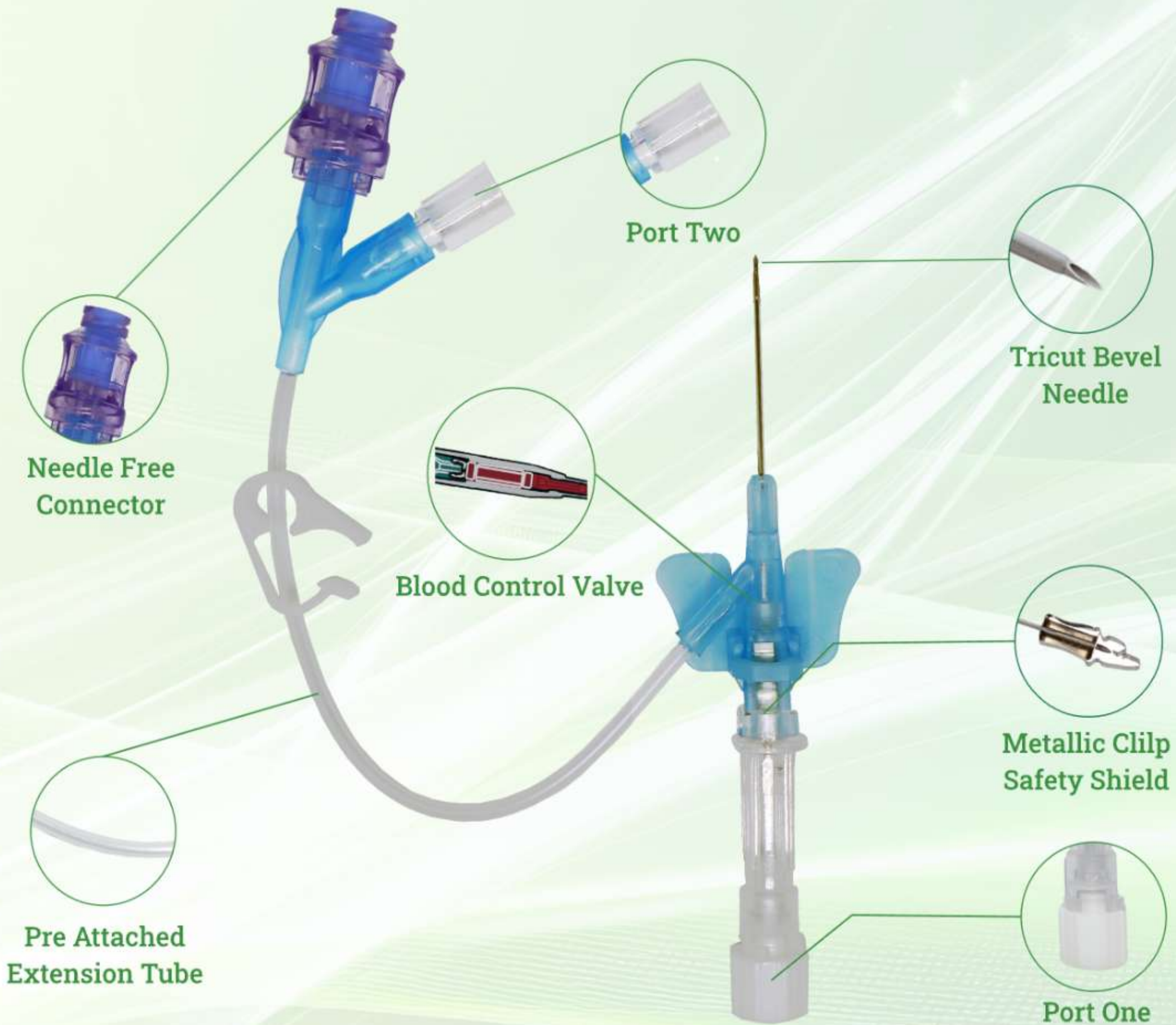

Technical Specification

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water Flow (ml/min)
18	Green	0.9-1.3	45	90
20	Pink	0.8-1.3	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Gen-X	200	72.50	34.00	15.20

Dial Flow Regulator



Flow Regulator with Y-Site

- Designed to regulate the flow of IV fluid from an infusion set into a IV Cannula.
- Drum regulator adjustment from 5 ml/ hr ~ 250 ml/hr.
- Built- in latex free Y-injection Site.
- Compatible to various devices with male & female Luer Lock
- Made up of medical grade non toxic and non kink-able, soft PVC tubing.


Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Flow Regulator	200	38.00	29.00	42.00

Pro Flow Delta™

Flow Regulator with Needle Free Connector

- Prevent needle stick injury by eliminating needle for IV access.
- Reduce infection risk by reducing entry points for bacteria.
- Neutral or positive displacement connector are designed to prevent blood reflux into the catheter.
- Drum regulator adjustment from 5 ml/ hr ~ 250 ml/hr.
- Made up of medical grade non toxic and non kink-able, soft PVC tubing.

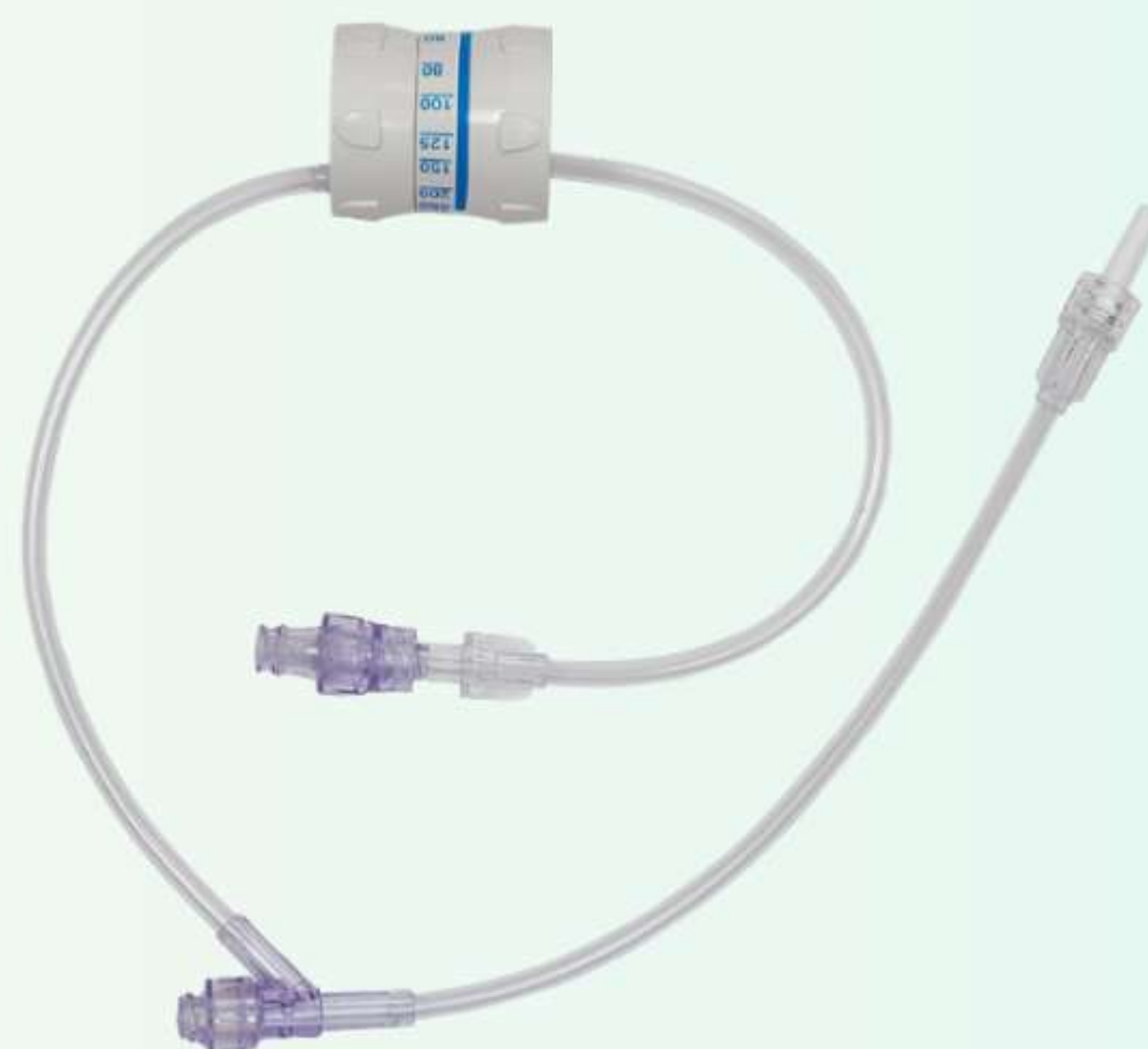

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Pro Flow Delta	200	38.00	29.00	42.00

3 WAY STOP COCK & EXTENSION

Proway

Three Way Stop Cock

- A Three Way Stop Cock facilitates multiple lines through single IV access.
- It is made up of a transparent medical grade polycarbonate.
- It supports smooth tap rotation at 360 degrees.
- Can be used with intra-arterial and intra-venous systems.
- Uniform infusion flow without any turbulence or air bubble.
- Accurate and continuous flow of channels.
- Available in lipid resistant/non-lipid resistant.
- Available in 3 bar and 5 bar pressure
- Available with 45° click tap movement.



360° Rotating Tap




Available in Red Tap



Available in White Tap



Available in Red & Blue Pegs

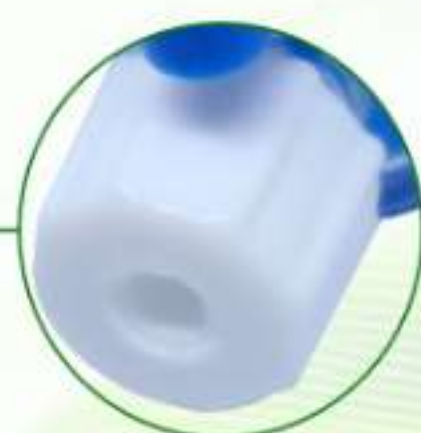

6% taper Luer Lock Connection

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Proway	500	54.00	29.00	18.00



Pro Ext Extension Tube

- Minimizes mechanical irritation of blood vessel being used for access.
- Reduces thrombotic events such as thrombophlebitis.
- Increase patient comfort.
- Ideal for administering contrast medium fluids during imaging procedure that require patient to enter a gantry e.g. CT & MRI scans.
- High quality Kink-Resistant tubing with smooth inner surface for laminar, non turbulent flow.
- Available in various tube lengths based on customer requirement.

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Pro Ext	200	38.00	29.00	42.00

Proway - Ext

3 Way Stop Cock with Extension Tube

- It combines the functionality of a three-way stopcock and an extension line.
- Minimizes mechanical irritation of blood vessel being used for access.
- Reduces thrombotic events such as thrombophlebitis.
- Increase patient comfort.
- Ideal for administering contrast medium fluids during imaging procedure that require patient to enter a gantry e.g. CT & MRI scans.
- High quality Kink-Resistant tubing with smooth inner surface for laminar, non turbulent flow.
- Available in various tube lengths based on customer requirement.



Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Proway Ext	200	38.00	29.00	42.00

AVAILABLE IN DEHP FREE

INFUSION SET

Pro-Inject

Built-in Air Vent, Bulb Latex

- Available with or without an air vent for compatibility with glass or plastic bottles.
- Transparent, flexible drip chamber for easy visualization of fluid flow and drop count.
- Made up of medical grade non toxic PVC.
- Roller Clamp facilitates precise flow rate control of IV fluids.
- 6% Universal Luer Lock/Slip fitting ensures leak-proof connection with IV cannulas.
- Medical-grade natural rubber latex (or latex-free version on request) for safe drug administration.
- Sharp, non-vented or vented spike for easy and safe penetration of IV fluid bottles.
- Sterilized by ethylene oxide (ETO) and individually packed to ensure patient safety.
- Customization available.

AVAILABLE IN DEHP FREE



Kink-Resistant Tube



Air Vent With Collar/Wings



Bulb Latex



Smooth Flow Controller



Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Pro Inject	800	66.50	38.00	37.00



Pro-Venex

Non Vented, Tube Latex

- Transparent, flexible Non-vented drip chamber.
- Made up of medical grade non toxic PVC.
- 6% Universal Luer Slip fitting ensures leak-proof connection with IV cannulas.
- Medical-grade natural rubber latex.
- Sharp, non-vented spike for easy penetration.
- Sterilized by ethylene oxide (ETO).
- Customization available.

IV Set For Oncology

- Cytotoxic & chemotherapy Drug Compatibility:
- Non-Toxic, DEHP-Free PVC Tubing
- Built-in Air Vent Spike
- Precision Flow Regulator:
- 6% Universal Luer C connector.
- Latex-Free Components.
- Sterile & Pyrogen-Free.



IV Set with Y-Site

- Y-Injection Site allows the administration of additional drugs during IV infusion without disconnecting the line.
- Built-in Air Vent Spike
- Precision Flow Regulator:
- 6% Universal Luer connector.
- Latex-Free Components.
- Sterile & Pyrogen-Free.

AVAILABLE IN DEHP FREE

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Pro-Venex	1000	66.50	38.00	37.00
IV Set - Oncology	500	66.50	38.00	37.00
IV Set with Y Site	500	66.50	38.00	37.00

SAFETY INFUSION SET

Pro-Inject Safe

Built-in Air Vent, Y-Site, Rotating Luer Lock and Priming Cap

- The Priming Cap ensures Aseptic Priming and enhances patient safety.
- Rotating Male Luer Lock ensures a Secure and Leak-Proof connection during infusion.
- Y-Injection Site allows the administration of additional drugs during IV infusion without disconnecting the line.
- Built-in Air Vent Spike
- Precision Flow Regulator:
- 6% Universal Luer connector.
- Latex-Free Components.
- DEHP Free.
- Sterile & Pyrogen-Free.



Packaging Details

Product	No. of Pieces Carton	Dimension (LxWxH in cm)		
		L	W	H
Pro Inject Safe	400	66.50	38.00	37.00

IV Set with Air Stop Filter



- 3-part drip chamber design with Air Stop Filter prevents air embolism, ensuring that air bubbles do not reach the patient during infusion.
- Roller Clamp enables precise control of fluid flow during infusion.
- 6% Luer Connector ensures secure & leak-proof connection.
- Made from non-toxic, DEHP-free, and biocompatible materials suitable for all types of IV fluids.
- Reduces risk of complications during high-risk or long-duration IV therapy.
- Sterile and Pyrogen-Free.

IV Set with Needle Free Y-Site


- Needle-Free Y-Site Port enables drug or fluid administration without puncturing the port with a needle, ensuring needle-stick injury prevention and aseptic drug delivery.
- Maximum safety and convenience.
- Built-in Air Vent Spike
- Precision Flow Regulator:
- 6% Universal Luer connector.
- Latex-Free Components.
- DEHP Free.
- Sterile & Pyrogen-Free.



Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
IV Set Air Stop Filter	400	66.50	38.00	37.00
IV Set with Needle Free Y-Site	400	66.50	38.00	37.00

Blood Transfusion Set

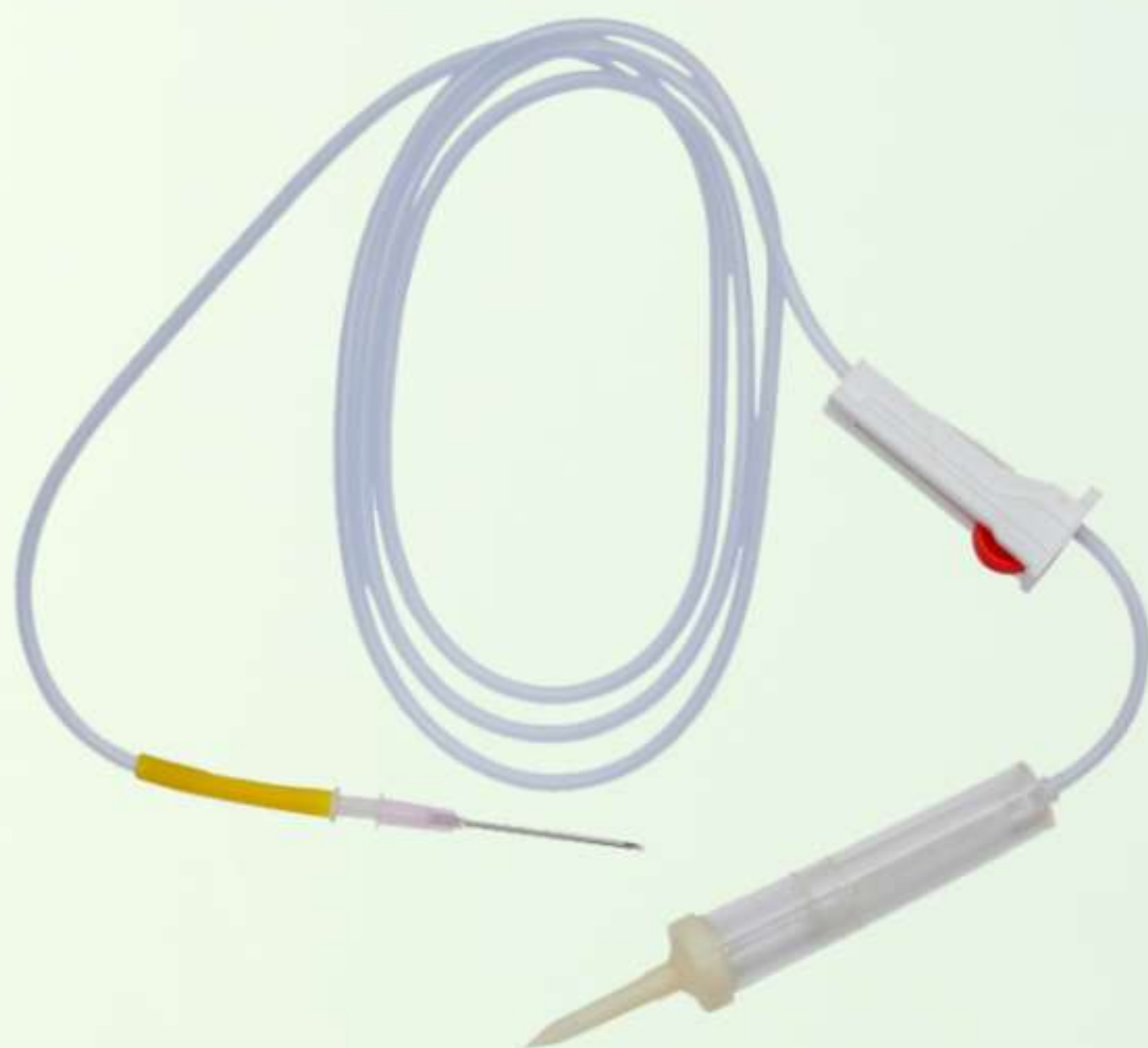

Pro-Fusion

Non Vented

- Non-Vented Spike is designed specifically for use with collapsible blood bags.
- Eliminate the need for external air venting and reduce contamination risks.
- Drip Chamber with Filter to trap clots and particulate matter.
- The Transparent Drip Chamber provides clear visibility of blood flow and drop rate.
- Flexible, kink-resistant tube made from medical-grade PVC to ensure compatibility with blood components.
- Roller Clamp facilitates precise control of flow rate during transfusion.
- 6% Luer Connector ensures a leak-proof and secure connection.
- Sterile & Pyrogen-Free.

Pro-Fusion Delta

Built-in Air Vent

- Equipped with a built-in air vent and hydrophobic filter that allows smooth flow of blood from glass or rigid bottles, while preventing fluid backflow or contamination.
- Drip Chamber with Filter to trap clots and particulate matter.
- The Transparent Drip Chamber provides clear visibility of blood flow and drop rate.
- Flexible, kink-resistant tube made from medical-grade PVC to ensure compatibility with blood components.
- Roller Clamp facilitates precise control of flow rate during transfusion.
- 6% Luer Connector ensures a leak-proof and secure connection.
- Sterile & Pyrogen-Free.




AVAILABLE IN DEHP FREE

Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Pro Fusion	400	66.50	38.00	37.00
Pro Fusion Delta	400	66.50	38.00	37.00

Measured Volume Set

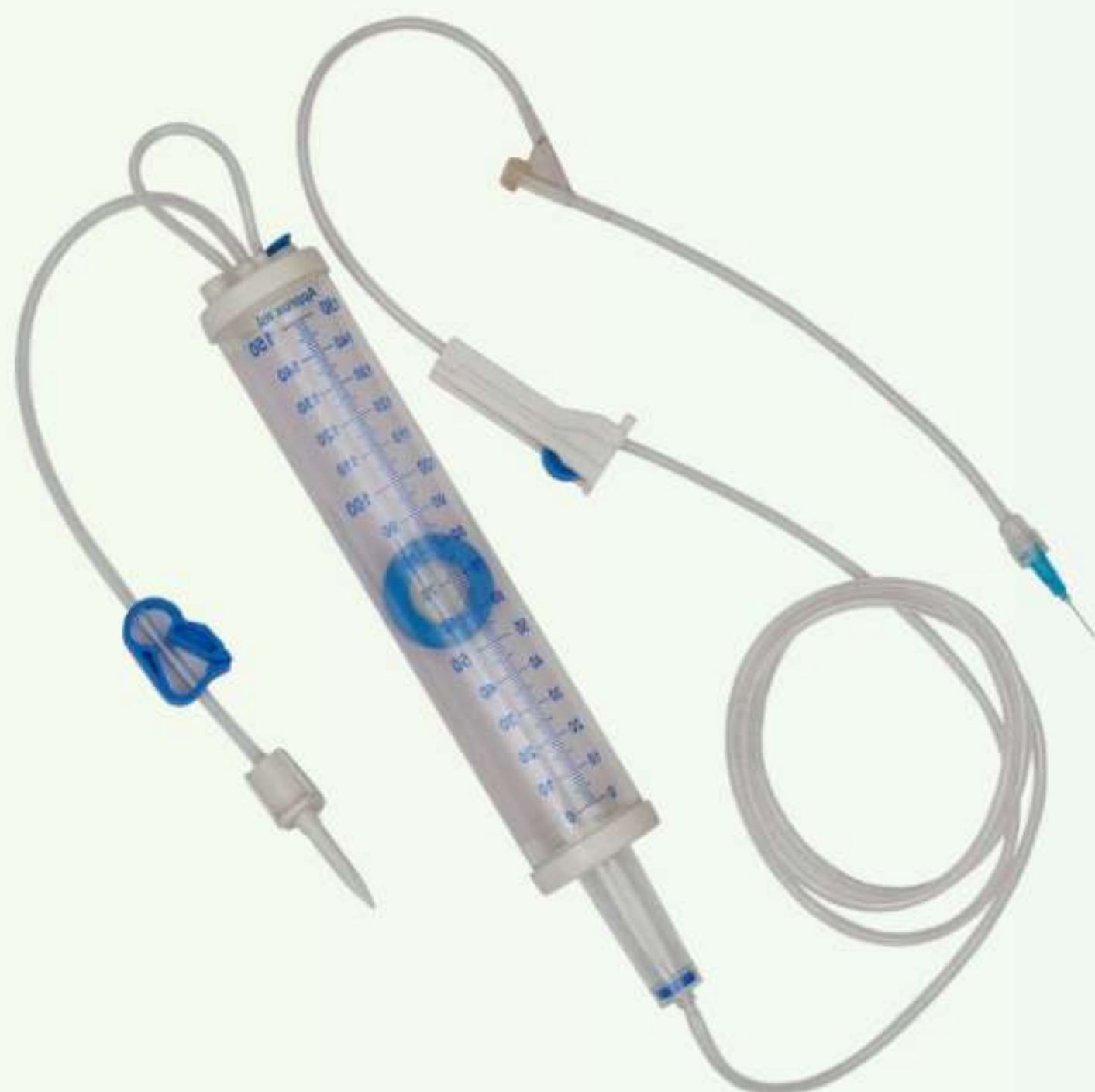

For Fluid

- Built-in Air Vent chamber for precise measurement and delivery of IV fluids.
- Y-Injection Port allows for the administration of additional medication.
- Vented spike allows air to enter the IV fluid container, ensuring a smooth and uninterrupted flow of fluid during infusion.
- Automatic shut-off floating valve prevents air trapping and also helps to indicate the volume of fluid in the chamber.
- Roller Clamp provides accurate control of flow rate.
- Soft & Kink-Resistant Tubing made from medical-grade, non-toxic PVC.
- 6% Luer Connector ensures a secure and leak-proof connection.
- Delivers approximately a flow rate of 60 drops/ml.
- Sterile & Pyrogen-Free.

AVAILABLE IN DEHP FREE

For Oncology

- DEHP-free Amber-Colored Tubing and Chamber protect medications from UV and visible light exposure, maintaining their potency during administration.
- Suitable for oncology and critical care drug infusion.
- Y-Injection Port allows for the administration of additional medication.
- Vented spike allows air to enter the IV fluid container, ensuring a smooth and uninterrupted flow of fluid during infusion.
- Automatic shut-off floating valve prevents air trapping and also helps to indicate the volume of fluid in the chamber.
- Sterile & Pyrogen-Free.



Packaging Details

Product	Size	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
			L	W	H
MV Set	150 ML	100	53.50	32.00	42.50

Hypodermic Syringe



Luer Lock

- Luer Lock Connector provides a twist-and-lock secure connection between syringe and needle, preventing accidental detachment.
- High-Precision Graduations with clear, accurate markings on the syringe barrel.
- Made from transparent, non-toxic medical-grade material.
- Smooth Plunger Action for effortless and controlled fluid aspiration or injection.
- Compatible with all standard Luer Lock needles of various gauges and lengths.
- Sterile & Single-Use.

Packaging Details

Product	Size	No. of Pieces /Carton	Dimension (LxWxH in cm)		
			L	W	H
Syringe Luer Lock	2,3ml	2000	58.50	26.50	48.50
	5ml	2000	66.00	30.50	49.00
	10ml	1800	58.00	43.50	44.50

Luer Slip

- Luer Slip Connector provide an easy push-on connection for needles, allowing quick assembly and smooth operation.
- High-Precision Graduations with clear, accurate markings on the syringe barrel.
- Made from transparent, non-toxic medical-grade material.
- Smooth Plunger Action for effortless and controlled fluid aspiration or injection.
- Compatible with all standard Luer Lock needles of various gauges and lengths.
- Sterile & Single-Use.

Packaging Details

Product	Size	No. of Pieces /Carton	Dimension (LxWxH in cm)		
			L	W	H
Syringe Luer Slip	2,3ml	2400	56.50	46.00	33.50
	5ml	1800	56.50	46.00	33.50
	10ml	800	56.00	40.00	28.00



LATEX GLOVES



Hand Care

Surgical Gloves

- Made from high-quality, medical-grade natural latex, offering superior elasticity and strength.
- Available in Powder-Free for reduced allergenic reactions or lightly Powdered for ease of donning.
- Ergonomic design for comfortable fit, flexibility, and reduced hand fatigue.
- Smooth Fingertips ensures enhanced grip and precision when handling surgical instruments.
- High Barrier Protection against bloodborne pathogens, body fluids, and cross-contamination.
- Latex Elasticity provides a snug fit, excellent stretch, and superior tactile sensitivity for delicate procedures.
- Sterile & Single-Use.
- Available Sizes: 6.0, 6.5, 7.0, 7.5, 8.0, 8.5.

Packaging Details

Product	Size	No. of Pieces /Carton	Dimension (LxWxH in cm)		
			L	W	H
Surgical Gloves	6.0 - 8.5	500	66.50	25.00	44.50

Hand Care

Examination Gloves

- Made of medical-grade natural latex provides excellent elasticity, strength, and durability.
- Available in Powder-Free for reduced allergenic reactions or lightly Powdered for ease of donning.
- Textured Fingertips ensure enhanced grip for handling instruments, equipment, or patient interactions.
- Barrier Protection against blood, body fluids, and contaminants, reducing the risk of cross-infection.
- Ideal for examinations, diagnostics, and routine procedures.
- Single-Use, Disposable & non-sterile.
- Available Sizes: XS, S, M, L, XL.



Packaging Details

Product	Size	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
			L	W	H
Examination Gloves	XS to XL	1000	33.50	22.50	21.00

RESPIRATORY RANGE



Oxygen Mask

- Made of soft, flexible, and non-toxic medical-grade PVC material ensures comfort and safety for patients.
- Ergonomic Design fits snugly over the nose and mouth for efficient oxygen delivery without leakage.
- Adjustable Elastic Strap allows a secure and comfortable fit.
- Transparent Mask Body ensures visual monitoring of the patient's breathing and oxygen flow.
- Compatible with standard oxygen supply tubing for seamless connection.
- Single-Use & Disposable.
- Available Sizes: Adult, Pediatric, & Neonatal.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
			L	W	H
Oxygen Mask	All	200	57.00	37.00	38.00


Nebuliser Mask

- Made of soft, flexible, and non-toxic medical-grade PVC material ensures comfort and safety for patients.
- Adjustable mask design provides a secure and comfortable fit, reducing medication leakage.
- Efficient Aerosol Delivery ensures maximum inhalation of medication with minimal wastage.
- Compatible with Standard Nebulisers.
- Single-Use & Disposable.
- Available Sizes: Adult, Pediatric, & Neonatal.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
			L	W	H
Nebuliser Mask	All	200	55.00	48.00	46.00



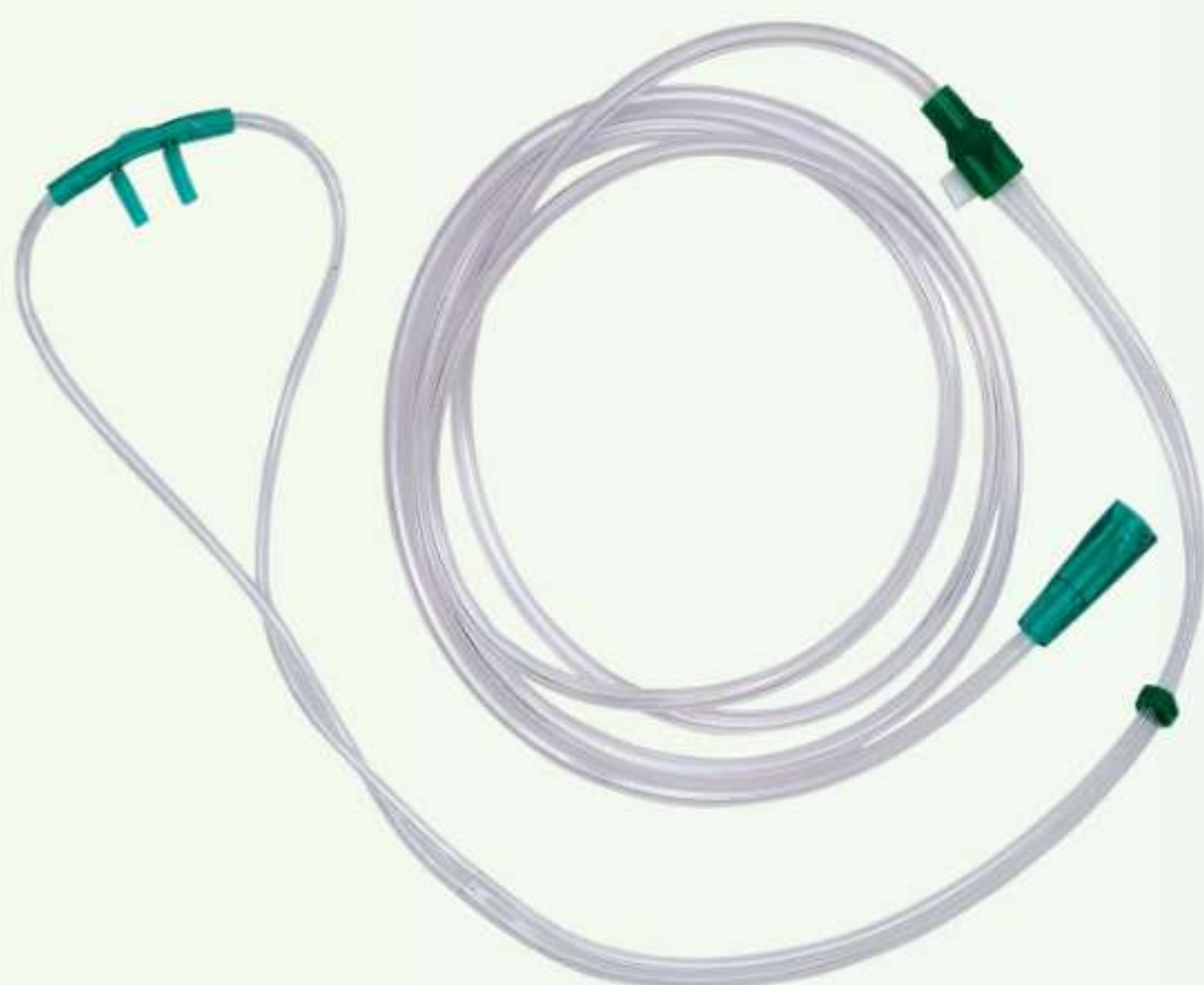

Twin Bore Nasal Catheter

- Dual-lumen (Twin-Bore) tubing delivers consistent oxygen flow and minimises resistance.
- Soft, Flexible Nasal Prongs ensure a comfortable fit in the patient's nostrils.
- Made of Non-toxic, soft medical-grade PVC ensures safety and durability.
- Compatible with standard oxygen flow meters for precise oxygen administration.
- Single-Use & Disposable.
- Available Sizes: Adult, Pediatric, & Neonatal.

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
			L	W	H
Twin Bore Nasal	All	200	45.00	40.00	48.00

Prociser

3 Ball Lung Exerciser

- Three-Ball System: Each ball represents a different resistance level, allowing users to measure inhalation strength and track progress over time.
- Transparent chamber allows users to monitor each ball's rise, providing immediate feedback on breathing performance.
- Made from non-toxic, lightweight & durable Medical-Grade plastic, ensuring safety, hygiene, and long-term use.
- Supports gradual respiratory training, with increasing resistance as lung strength improves.
- Portable and Easy to Use.
- For home use, hospital therapy, or rehabilitation programs.



Packaging Details

Product	No. of Pieces /Carton	Dimension (LxWxH in cm)		
		L	W	H
Prociser	28	58.00	37.00	38.00

OTHER PRODUCTS

Medulla

Spinal Needle

- High-quality stainless steel Cannula provides excellent strength and flexibility for accurate needle placement and reduced risk of breakage.
- Precisely Ground Quincke-Type Bevel Tip ensures sharp, smooth penetration with minimal tissue trauma and easy identification of dura puncture.
- Transparent Hub allows visual confirmation of cerebrospinal fluid (CSF) backflow, ensuring proper needle positioning.
- 6% Luer connector ensures secure connection with syringes and other compatible devices.
- Fitted snugly to prevent tissue coring and maintain needle patency during insertion.
- Size identification by Colour Coding facilitates quick and easy size selection during procedures.
- Sterile and Single-Use.
- Available Sizes: 15G, 16G, 17G, 18G, 19G, 20G, 21G, 22G, 23G, 24G, 25G, 26G, 27G.



Technical Details

Gauge	Colour	OD (in mm)	Effective Length (± in mm)
15	Blue-Gray	1.8	90
16	White	1.6	
17	Red-Purple	1.4	
18	Pink	1.2	
19	Cream	1.1	
20	Yellow	0.9	
21	Dark Green	0.8	
22	Black	0.7	
23	Deep Blue	0.6	
24	Purple	0.55	
25	Orange	0.5	
26	Brown	0.45	
27	Med-Gray	0.4	

Packaging Details

Product	Size	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
			L	W	H
Medulla	15G to 27G	500	45.00	36.00	18.00

Yankauer Suction Set

- Designed for effective removal of secretions, blood, and fluids from the surgical site or airway
- Transparent Yankauer handle is made from clear, rigid medical-grade plastic for easy visibility of fluids during suctioning.
- Smooth and rounded tip to prevent tissue trauma and ensure gentle suctioning.
- Kink-resistant tube made from non-toxic, medical-grade PVC ensures steady suction flow and user comfort.
- Universal Connector fits securely with standard suction apparatus or tubing systems.
- Sterile and Single-Use.



Injection Stopper/Heparin Cap

- Designed to maintain the sterility of IV lines, prevent leakage, and provide safe closure between infusions.
- Leak-proof seal provides a secure closure to prevent contamination, leakage, and air ingress in IV lines.
- Smooth, Soft Injection Membrane allows easy and repeated needle punctures without particle shedding or leakage.
- 6% Tapered Luer fitting.
- Sterile and Single-Use.

Combi Luer Lock

- Designed to fit all standard Luer Lock and Luer Slip connectors used in IV therapy and transfusion systems.
- Leak-proof design provides a secure, airtight seal to prevent leakage or air entry during fluid administration.
- Manufactured from transparent, non-toxic, and pyrogen-free medical-grade PVC, ensuring safety and durability.
- Available in Red & Blue Colours.
- Sterile and Single-Use.



Packaging Details

Product	No. of Pieces/ Carton	Dimension (LxWxH in cm)		
		L	W	H
Yankauer Suction	125	60.00	36.00	42.00
Injection Stopper	2500	54.00	29.00	18.00
Combi Luer	2500	54.00	29.00	18.00



Corporate Office

📍 136, 2nd Floor, Avtar Enclave
Paschim Vihar, New Delhi -110063, India
☎ +91-11-4573-3068
📱 +91-9910014706
✉ sales@larsmedicare.com

Americas Headquarters

📍 555 Industrial Dr. Milton,
ON Canada, L9T 8K1
📱 +1-905-202-0305
📱 +1-647-299-5115
✉ sameer@larsmedicare.com

Production Facilities

📍 Kila No. 16-17, Sultanpur
Opp. Sports Authority of India,
Near Bahalgarh Chowk,
Sonapat, Haryana - 131021, India

📍 Survey No: 801,
Chhadwada, Bhachau, Kutch,
Gujarat – 370140,
India

www.larsmedicare.com 🔍

- The description and specification in the brochure are correct at the time of printing.
- In pursuit of continual improvement, we reserve the right to change the design and specification of any of our products without prior notice.
- All photographs, product designs and trademarks represented in the brochure are the sole proprietary intellectual property of Lars Medicare Pvt. Ltd.

Catalogue Number:
LM/CAT/PD/2025/07